



## Installing OilPro V3.2.17 on Windows 8:


1. Insert the CD and browse to the folder Installation
2. The contents will be:

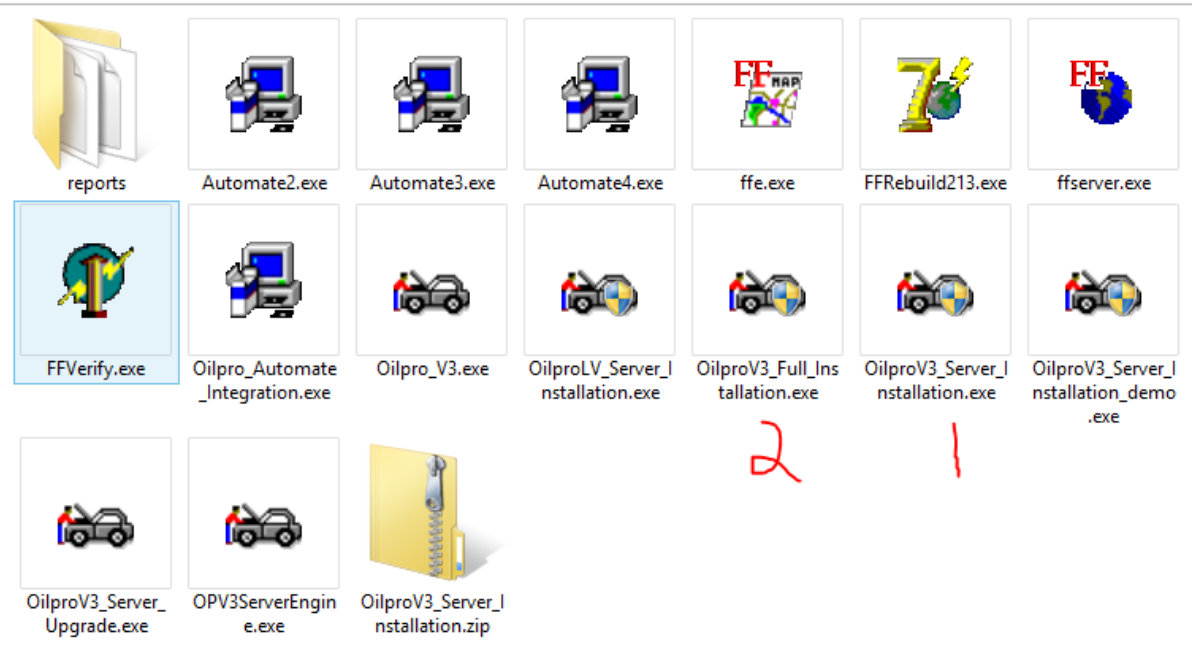

3. Right click on OilproV3\_Server\_Installation.exe and “Run as administrator”
4. Make sure the install path is to C:/Program Files/Oilpro\_Client

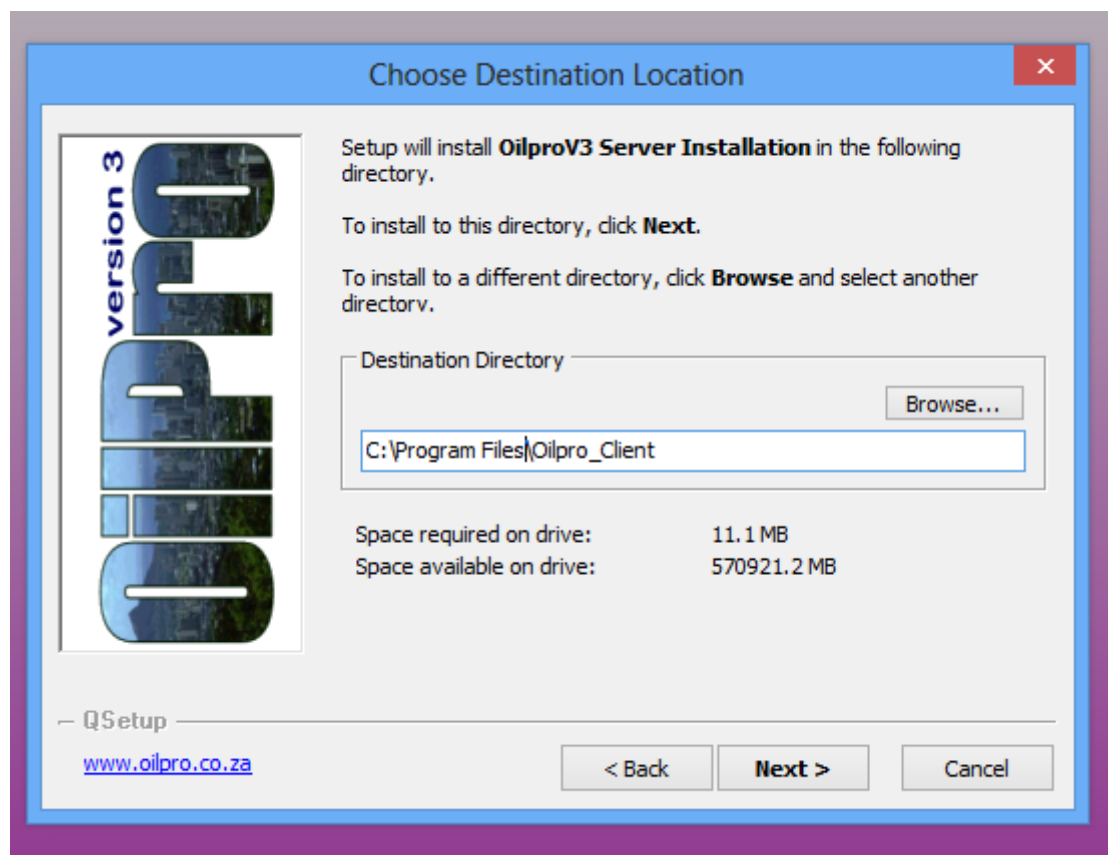


5. and not C:/Program Files (x86)/Oilpro\_Client (simply delete the (x86) so it looks like above.

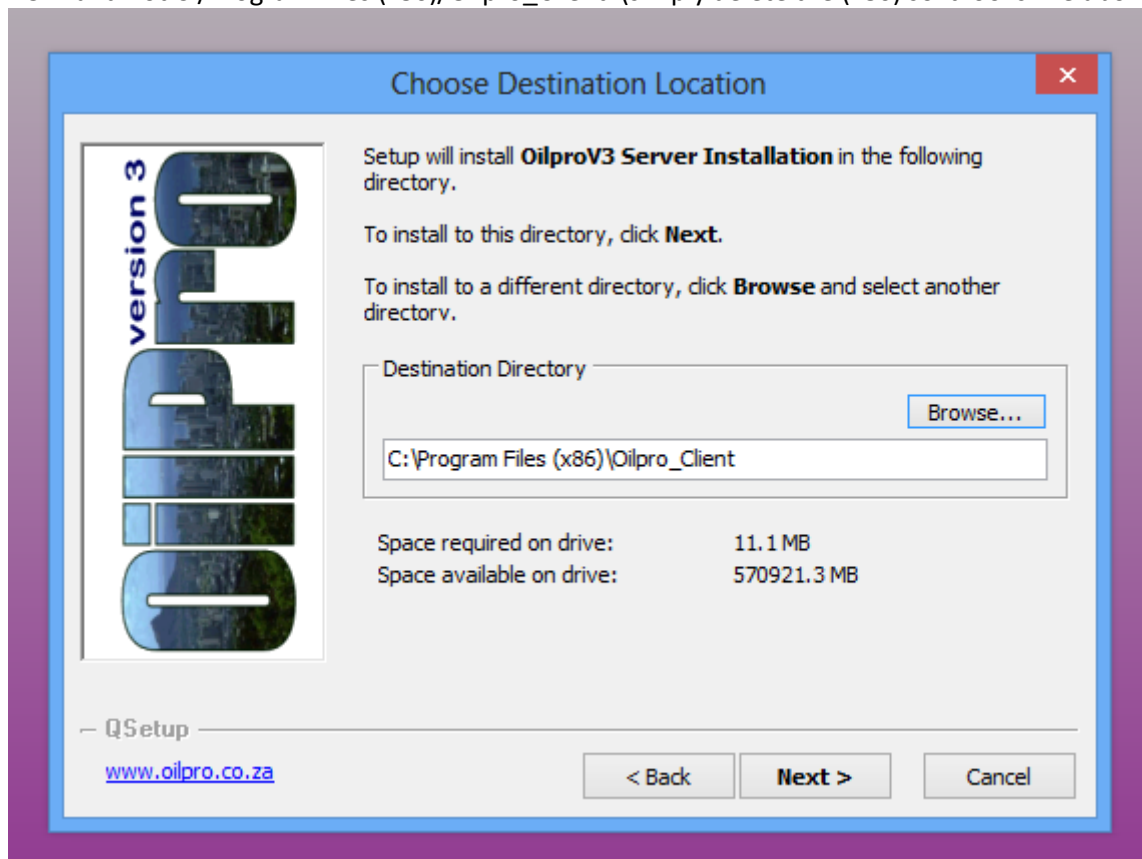

6. Do the same with the OilPro Client (OilproV3\_Full\_Installation.exe) and make sure you install it by right clicking and selecting "Run as Administrator"
7. Once both are installed browse to C:/Program Files (x86) and cut the Oilpro\_Data folder out of there and paste it to C:/Program Files.
8. Both Oilpro\_Data and Oilpro\_Client MUST reside under C:/Program Files.

## Installing the device driver on Windows 8:

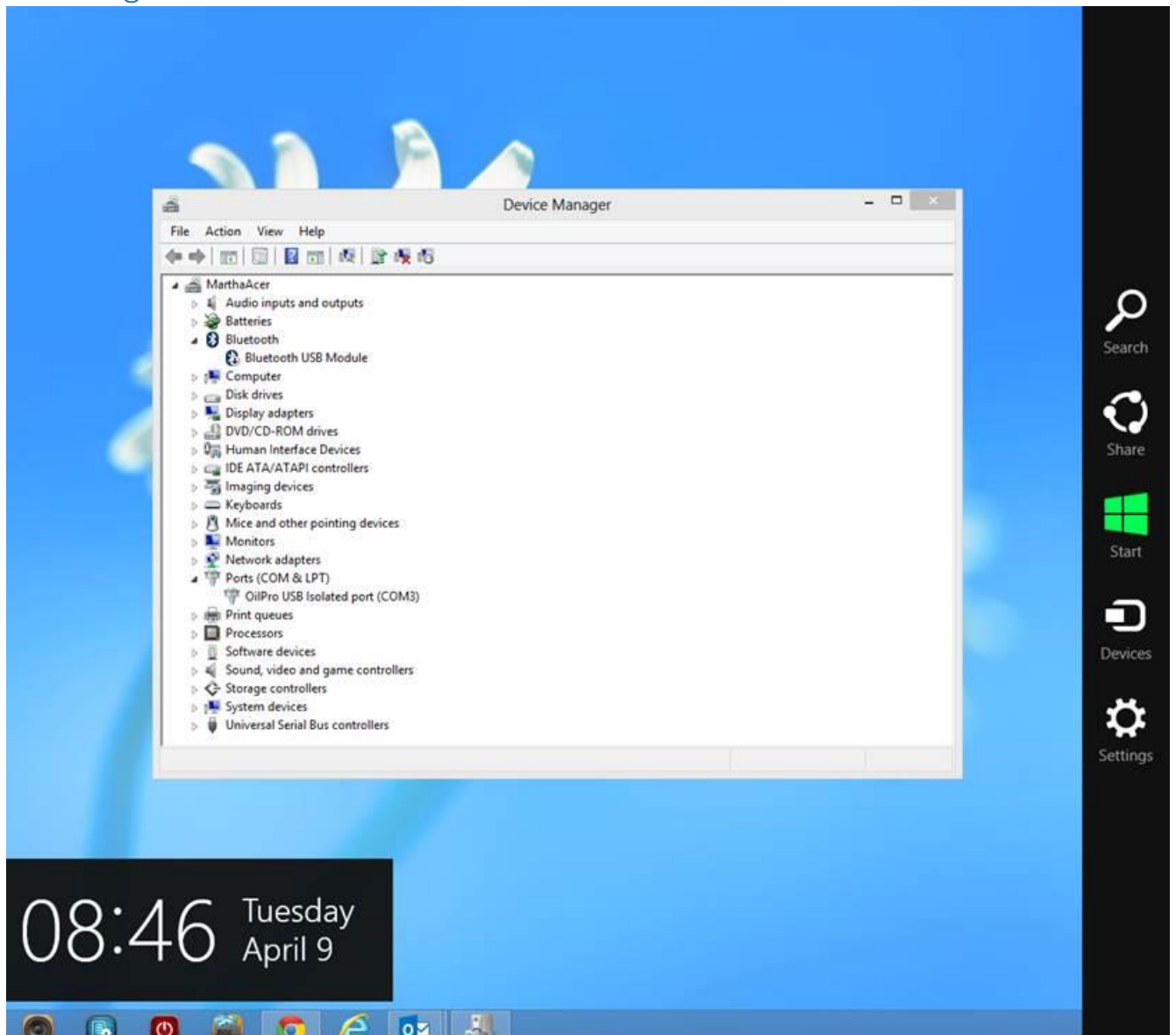

I managed to install it by following these steps:

Disable driver signature enforcement in Windows 8 (Works with 64-bit)

Select Settings (right bar)/More PC Settings (at bottom)

Then go General/Advanced Startup/Restart Now/Troubleshoot/Advanced Options/Windows Startup Settings/Restart

As your PC is starting up be pressing F8 to get to the menu where you can choose to Disable driver signature enforcement.

Remember to make sure the port number assigned by Windows matches the port number under the OilPro Server Local settings.